

ENDO 3

Ultrasonic Endo Activate Device



Versatile | User-friendly | Smart



User-friendly ring switch button

Three power levels

Three power levels available:
E1, E2, E3, power from low to high in order

Four time modes

Four time modes optional:
"10, 15, 20, --"

The working time in the irrigation mode:
"10s, 15s, 20s, 10min"

Highly convenient design More efficient irrigation

Smart microcomputer control

Microcomputer automatic control is adopted. The high-frequency oscillation signal is generated by the high-frequency oscillation circuit and acts on the handpiece head, making the tip resonate by the ultrasonic vibration stimulation.

Wireless control handy to operate

Magnetic suction contact charging
Wireless operation
Compact and light, easy for irrigation



Clear LCD screen easier observation

Visually display the machine
working state



User-friendly ring switch button

Ring switch button to control the
machine

Facilitate the treatment of different
angles and different tooth positions



Stable performance

Appropriate vibration amplitude

Tips with low vibration and high efficiency
More efficient and reliable irrigation




Stable performance Intelligent tracking

Fully automatic frequency tracking system
Automatically search for best working state



Compatible tips:

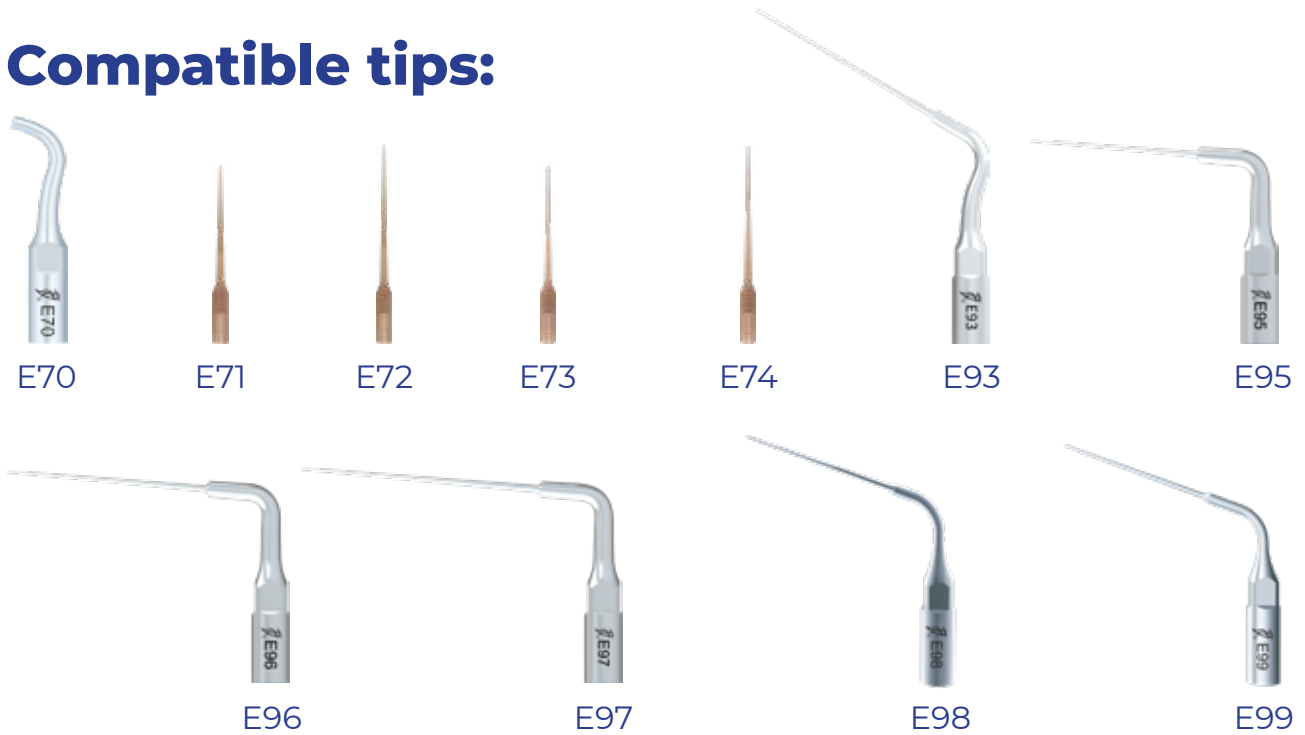

Standard configuration: E96×2, E99×2, E70, E71×10

High configuration: E96×2, E99×2, E70×2, E71×10, E72×10

Compatible tips: E70, E71, E72, E73, E74, E93, E95, E96, E97, E98, E99

To canal irrigation after the root canal preparation

E70: Endochuck, used for install E71, E72, E73 o E74 r root canal cleaning or irrigation

E71, E72, E73, E74: Polymer material without cutting ability

E93: 20#, taper 1%, tip end length 21.5mm, Ni-Ti

E95: 20#, taper 2%, tip end length 18mm, stainless steel

E96: 20#, taper 2%, tip end length 21mm, stainless steel

E97: 20#, taper 2%, tip end length 25mm, stainless steel

E98: 25#, taper 2%, tip end length 18.5mm, titanium alloy

E99: 25#, taper 2%, tip end length 18.5mm, titanium alloy



Guilin Woodpecker Medical Instrument Co., Ltd.
Information Industrial Park, Guilin National Hight-Tech
Zone, Guilin Guangxi, 541004 P. R. China

Wholesale importer Woodpecker / DTE
for Spain and Portugal

RUMAR Cedeira S.L.

C. del Obispo Golfín, 11, Bajo A, 28430 Alpedrete, Madrid

Phone.: 91 849 01 04

Cell Phone / Whatsapp: +34 626 327 863

Always in stock! Delivery in 24 hours.

Devices, accessories, spare parts and after-sales service

Your trusted distributor:

RUMAR
mayorista dental