



The choice of millions
of dentists worldwide

U600

Ultrasonic Scaler

Pleasant & Efficient



Intelligent power output
Experience painless treatment

U600 Ultrasonic Scaler

A global professional scaler manufacturer focuses on technological innovation for 30 years.

International high-tech enterprise
30 years innovation
Woodpecker scaler owns 101 national patents
9 national invention patents

Effective and thorough endodontic irrigation



After irrigation, there is no smear layer on surface but a few amount of dentin debris in the root canal.

Root canal irrigation is the most widely used clinical support methods. Its purpose is to remove the substance inside root canal, dissolve diseased tissue, remove pathogenic microorganisms, and remove smear layer. Lubricate the root canal internal wall. During root canal preparation process, prevent the infectious substance from being pushed to the deep part or out of the apical foramen, reducing the risk of instrument fracture in root canal.


Compared with the traditional syringe irrigation, ultrasonic irrigation can better remove bacteria, dentin debris, and smear layer in root canal. And it can well clean special structures such as the lateral branch canal and root canal narrow area. Currently, ultrasonic irrigation is widely used in endodontic therapy in order to receive a better root canal preparation effect.



Sine wave drive

Avoid instantaneous stimulation caused by power surge, ensuring comfortable treatment.

U600 Mild power output

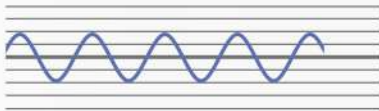

Rectangular power out of other brand scaler



High-quality handpiece

Aluminum shell can be autoclaved, not easy to cause cracking or deformation. Latest handpiece sealing technique enables long durability. User-friendly design brings better hand feeling.



Automatic water supply system

Easy to change water, eliminating redundant steps and increasing the working efficiency. Solutions such as hydrogen peroxide, sodium hypochlorite, and chlorhexidine can be used. Clear 50ml scale mark makes it convenient to formulate solutions.



> High configuration

G4	G12	P1	P3	P8	E60	E62
x2	x2	x2	x2	x2	x2	x2

Standard configuration

G4	G12	P1	P3	P8	E60	E62
x1	x1	x1	x1	x1	x1	x1

Note: Woodpecker reserves the right of modifying the tip configuration without further notice.


Supragingival scaling



Implant maintenance



Subgingival scaling

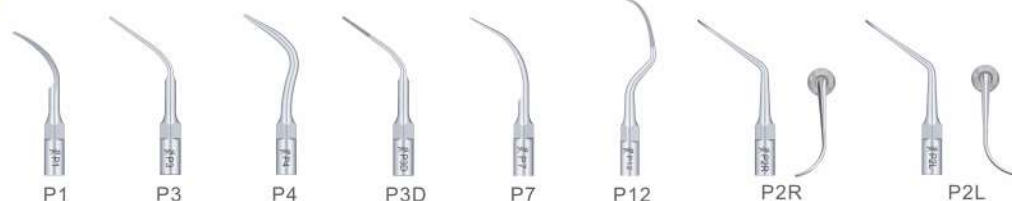


Root canal treatment

Endodontic irrigation



For removal of calcification and bad fillings in root canal

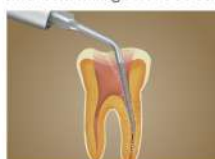

For removal of calcification at root canal orifice and top 1/3 of root canal



Root canal retreatment Removal of foreign body



Special material (E6, E7)
Thin, flexible, good
conductivity to ultrasound

For removal of
foreign body
at bottom 1/3
of root canal

For removal of
foreign body
at top 1/3 of
root canal

For removal of
foreign body and
instrument fractures
(stainless steel)
in root canal

For removal of
temporary cement and
instrument fracture
(stainless steel)
in root canal

Locate the root canal orifice



Diamond-coated ball head
Locate root canal orifice with
security, accuracy and efficiency

Cavity preparation



G30, For preparation
between pits and fissures



G31, For the preparation
of mesioocclusion cavity



G32, For the preparation
of distal cavity



For preparation
between pits
and fissures



For the
preparation of
mesioocclusion
cavity



For the
preparation
of distal cavity



For preparation in
occlusal surface
and cervix dentis



For the
preparation of
mesial surface



For the
preparation of
distal surface



Guilin Woodpecker Medical Instrument Co., Ltd.
Information Industrial Park, Guilin National High-Tech
Zone, Guilin Guangxi, 541004 P. R. China

Importador mayorista Woodpecker / DTE
para España y Portugal

RUMAR Cedeira S.L.

C. del Obispo Golfín, 11, Bajo A, 28430 Alpedrete, Madrid

Tlfno.: 91 849 01 04

Tlfmo. móvil / Whatsapp: +34 626 327 863

¡Siempre en Stock! Entrega en 24 horas.

Dispositivos, accesorios, repuestos y servicio postventa

Su distribuidor de confianza:



RUMAR
mayorista dental